



## INTELECT® **MOBILE 2**

THE NEXT GENERATION  
IN ELECTROTHERAPY





## INTELECT® MOBILE 2

### THE NEXT GENERATION IN ELECTROTHERAPY

The Intelect® Mobile 2 is the next generation in physical medicine modalities, its intuitive design is clever and its features, usability and simplicity make Intelect Mobile 2 the choice for today's modern practice.

The Intelect Mobile 2 is simple to use and set up thanks to the built in suggested protocol settings developed around the current clinical practices.

The Intelect Mobile 2 has an easily accessible high resolution anatomical library which allows you to explain visually to the patient their condition and the treatment pathway.

Smart Bluetooth module and the intellect Connect App give the Intelect Mobile 2 device the ability to download data, produce reports and upload treatment data to additional devices within the practice. The Intelect Connect App will allow patient data to be collected and share between all Intelect Connect ready devices.





# KEY FEATURES

## INNOVATION AND SIMPLICITY THROUGH DESIGN



### LARGE GRAPHIC COLOR TOUCH SCREEN

High resolution 7" capacitive touch screen for access to intuitive user interface and extensive anatomical library. Suggestions for the placement of electrodes and areas for ultrasound applicator placement are in easy to use and find guidelines. The touch screen also gives real time information as to the progress of the individual treatment and a summary after treatment ends.



### NEW ULTRASOUND APPLICATOR

Ergonomic design fit for the hand Improved coupling detection and monitoring.

### ERGONOMIC DESIGN

New contemporary modern design to fit in any therapy office, but versatile enough with the integral option of the battery to make the device simple and convenient to use in the patients own home.



### MOBILE

Lightweight device with ergonomically designed handle and optional integral battery, eliminating the need to carry mains wires on home visits. Optional custom carry case available for transporting the device safely and securely.

### HOME USE

Intellect® Mobile 2 was designed to be compliant with IEC 60601-1-11 and can be used for treatment in the patients own home as well as the therapist office, giving flexibility and versatility for both clinician and patient alike.

### WINDOWS 10 APP

The addition of a smart Bluetooth module to the Intellect Mobile 2 and provision of the Intellect Connect App\* gives great flexibility in the handling and processing patient data, allowing the export of patient treatment data and protocols from device to PC allowing not only the ability to produce reports but also to share treatments and patient data between multiple Intellect Mobile 2 devices within a practice and in the near future with other smart enabled Chattanooga devices.



### OPTIONAL CART (OPTIONAL VACUUM UNIT)

Optional cart with storage drawers ensuring all accessories are convenient and easy to reach from the device, large wheels allow the cart and device to be moved easily around the office, add in the option of the vacuum unit to allow the use of vacuum in electrotherapy.



Download the Intellect Connect App free from [www.chattanooga rehab.com](http://www.chattanooga rehab.com)

# KEY FEATURES - INTERFACE

Intelect® Mobile 2 offers the next generation of intuitive and intelligent user interface - designed with the therapist in mind.

## INTELLIGENT INTERFACE

Buttons on the device and connection points for Modality accessories have interactive lights to show the status of the device or port.



## SUGGESTED PROTOCOL SETTING

Assists the user in setting the functions of the device and each therapy, settings are based on current clinical practice.

Suggested protocol settings (SPS) can be modified to produce custom treatment regimens and saved as a custom protocol.



## HI RESOLUTION ANATOMICAL LIBRARY

Illustrates an array of pathologies allowing simple explanation of conditions and treatment to the patient.



## FUTURE PROOF TECHNOLOGY

Intelect Mobile 2 has been designed with the ability of software upgrades allowing additional functionality to be added. The free software upgrades will be made available on the [chattanoogarehab.com](http://chattanooga rehab.com) website along with simple instructions on how to perform any upgrades and access to updated product manuals incorporating any software enhancements to the device.



# INDICATIONS

## STIM

### TENS

- SYMPTOMATIC RELIEF OF CHRONIC PAIN (1)
- MANAGEMENT OF POST-OPERATIVE PAIN (2)

### NMES

- MUSCLE RE-EDUCATION (3)
- INCREASING LOCAL BLOOD SUPPLY (4)
- RELAXATION OF MUSCLE SPASMS (5)
- MAINTAINING/INCREASING RANGE OF MOTION (6)

## ULTRASOUND

- RELIEF OF PAIN FROM MUSCLE SPASM (7)
- RELIEF OF PAIN FROM JOINT CONTRACTURE (8)
- RELIEF OF PAIN ASSOCIATED WITH LIGAMENT (9)
- SPRAINS, TENDINITIS AND MUSCLE SPASM (10)



## CLINICAL INFORMATION

1. The effectiveness of transcutaneous electrical nerve stimulation in the management of patients with complex regional pain syndrome: A randomized, double-blinded, placebo-controlled prospective study. Bilgili A, Çakır T, Doğan ŞK, Erçalık T, Filiz MB, Toraman F Nov 2016
2. Comparison of the Effect of Sensory-Level and Conventional Motor-Level Neuromuscular Electrical Stimulations on Quadriceps Strength After Total Knee Arthroplasty: A Prospective Randomized Single-Blind Trial. Yoshida Y, Ikuno K, Shomoto K.Dec 2017
3. Neuromuscular Electrical Stimulation Superimposed on Movement Early after ACL Surgery. Labanca L, Rocchi JE, Laudani L, Guitaldi R, Virgulti A, Mariani PP, Macaluso A. Mar 2018
4. A randomised controlled trial of electrostimulation effects on effusion, swelling and pain recovery after anterior cruciate ligament reconstruction: a pilot study. Ediz L, Ceylan MF, Turktas U, Yanmis I, Hiz O Oct 2011
5. A promising approach to effectively reduce cramp susceptibility in human muscles: a randomized, controlled clinical trial. Behringer M, Moser M, McCourt M, Montag J, Mester J. Apr 2014
6. Effect of electrical stimulation to long head of biceps in reducing gleno humeral subluxation after stroke. Manigandan JB, Ganesh GS, Pattnaik M, Mohanty P.Jan 2014
7. Thermal ultrasound decreases tissue stiffness of trigger points in upper trapezius muscles. Draper DO, Mahaffey C, Kaiser D, Eggett D, Jarmin J. Apr 2010
8. Overview of Therapeutic Ultrasound Applications and safety considerations. Miller,D., Smith, N., Bailey, M., Czarnota, G., Hynynen, k., & Makin I 2012
9. A comparison of the effectiveness of low-, moderate and high dose ultrasound therapy applied in the treatment of myofascial pain syndrome. Modern Rheumatology 24, 2014.
10. Efficacy of Therapeutic Ultrasound in Pain and Joint Mobility in Whiplash Traumatic Acute and Subacute Phases. Ruiz-Molinero C, Jimenez-Rejano JJ, Chillon-Martinez R, Suarez-Serrano C, Rebollo-Roldan J, Perez-Cabezas V Sep 2014




# INTELECT® MOBILE 2 RANGE



## STIM

The Intellect® Mobile 2 **STIM** is an electrotherapy system used with or without an optional cart, allowing for the inclusion of a Vacuum module. The Intellect Mobile 2 stim has available a selection of widely used waveforms allowing the therapist multiple treatment options, these waveforms are easily selected from the user interface and guidelines as to electrode placement are also given for the selected treatment area. This equipment is to be used only under the prescription and supervision of a licensed medical practitioner.

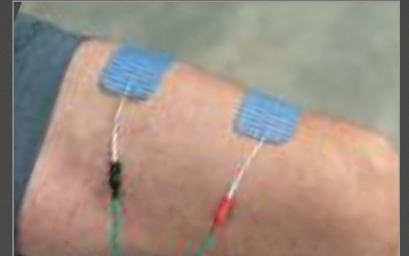

## ULTRASOUND

The Intellect Mobile 2 **ULTRASOUND** is an ultrasound therapy system used with or without an optional Cart. The Intellect Mobile 2 Ultrasound offers dual frequency Ultrasound at 1 or 3 MHz, Pulsed or continuous duty cycles (10%, 20%, 50% & 100%). Supplied with a 5cm<sup>2</sup> ultrasound applicator with additional optional heads available. Guidelines for treatment area are available on the user interface. Coupling detection can be monitored either audibly (set to on as default) or visually by the LED on the back of the applicator, treatment time indicator will stop when head is uncoupled so the of amount of ultrasound treatment delivered during the time period is accurately measured and reproducible. This equipment is to be used only under the prescription and supervision of a licensed medical practitioner.



## COMBO

The Intellect Mobile 2 **COMBO** is a two-channel electrotherapy, ultrasound therapy and Combo system used with or without an optional Cart, allowing for the inclusion of a Vacuum module. For combo therapy 5 waveforms are included for use: Interferential, VMS, High Volt, Asymmetric biphasic TENS and Symmetric biphasic TENS there is fully functional 1 and 3 MHz ultrasound frequencies and two channels of electrical stimulation output. This equipment is to be used only under the prescription and supervision of a licensed medical practitioner.



# TECHNICAL INFORMATION

INPUT: 100-240V AC, 1.0 to 0.42A, 50/60 Hz  
 ELECTRICAL CLASS: Class II  
 MODE OF OPERATION: Continuous



## ELECTRICAL TYPE (DEGREE OF PROTECTION)

ULTRASOUND: .TYPE BF  
 ELECTROTHERAPY: .TYPE BF  
 ELECTROTHERAPY VACUUM: .TYPE BF  
 ULTRASOUND & ELECTROTHERAPY: .TYPE BF

## ULTRASOUND SPECIFICATIONS

FREQUENCY: 1 MHz; 3 MHz  
 DUTY CYCLES: 10%, 20%, 50%, Continuous  
 PULSE REPETITION RATE: 16, 48, or 100 Hz  
 PULSE DURATION: 1 -31.25 ms  
 - Max (ON): 31.25 ms  
 - Min (OFF): 5 ms

## VACUUM SPECIFICATIONS

**POWER**  
 INPUT: 20-25 Vdc, maximum peak current 4A  
 ELECTRICAL TYPE: TYPE BF

## GENERAL CHARACTERISTICS

VACUUM RANGE: .0 to 600 mbar maximum +/-5%  
 VACUUM MODES: Continuous or Pulsed Continuous  
 10 setting over vacuum range, 60 mbar per setting,  
 +10 mbar to 10 mbar per setting



| SYSTEM DIMENSIONS & WEIGHT     | WIDTH      | DEPTH      | HEIGHT | WEIGHT (NO BATTERY) |
|--------------------------------|------------|------------|--------|---------------------|
| INTELECT MOBILE 2 HEAD UNIT    |            |            |        |                     |
| COMBO                          | 34cm       | 35.5cm     | 15cm   | 3.1kg               |
| ULTRASOUND                     | 34cm       | 35.5cm     | 15cm   | 2.8kg               |
| STIM                           | 25.5cm     | 35.5cm     | 15cm   | 2.9kg               |
| CART CONFIGURATIONS            |            |            |        |                     |
| Cart (Safe working load 6.5kg) | 48cm (max) | 52cm (max) | 96cm   | 10.1kg              |
| Cart with vacuum               | 48cm (max) | 52cm (max) | 96cm   | 11.5kg              |
| Device on cart                 | -          | -          | 111cm  | -                   |

| ULTRASOUND APPLICATOR FREQUENCY                      | 1 cm <sup>2</sup>  |                    | 2cm <sup>2</sup>   |                    | 5cm <sup>2</sup>   |                    |
|--|--------------------|--------------------|--------------------|--------------------|--------------------|--------------------|
|  | 1MHz               | 3MHz               | 1MHz               | 3MHz               | 1MHz               | 3MHz               |
| Effective Radiating Area ERA INTL (cm <sup>2</sup> ) | 1                  | 0.9                | 1.5                | 1                  | 2.5                | 2.7                |
| Max Output power in Continuous mode                  | 2W                 | 1.8W               | 3W                 | 2W                 | 5W                 | 5.4W               |
| Max Output power in Pulsed mode                      | 3W                 | 2.7W (*)           | 4.5W               | 3W                 | 7.5W               | 8.1W               |
| Max Amplitude in Continuous mode                     | 2W/cm <sup>2</sup> | 2W/cm <sup>2</sup> | 2W/cm <sup>2</sup> | 2W/cm <sup>2</sup> | 2W/cm <sup>2</sup> | 2W/cm <sup>2</sup> |
| Max Amplitude in Pulsed mode                         | 3W/cm <sup>2</sup> | 3W/cm <sup>2</sup> | 3W/cm <sup>2</sup> | 3W/cm <sup>2</sup> | 3W/cm <sup>2</sup> | 3W/cm <sup>2</sup> |

|                   |  |
|-------------------|--|
| Pulsed Biphasic   | VMS,VMS Burst,VMS FR,TENS Asymmetrical,TENS Symmetrical,HAN  |
| Pulsed Monophasic | High Volt,Rectangular,Triangular,Surged Rectangular,Surged Triangular, Micro Current, Trabert,Galvanic interrupted                     |
| Direct            | Galvanic, Low Level Galvanic   |
| Alternating       | IFC 4Pole,IFC Premod, IFC Isoplanar Vector, Russian, Diadynamic MF, Diadynamic DF, Diadynamic CP, Diadynamic LP, Diadynamic Cpid/Cpiso |

Chattanooga also offer online training for their customers with modules available for their modalities, please contact your Chattanooga representative for more information and sign up details.

\*Conformity :This device complies with the applicable standards EN/IEC 60601-1, CAN/CSA-C22.2 No.6011, UL Std. No 60601-1"  
 Classification according to MDD : class Ila device 93/42 60601-1  
 Legal manufacturer: DJO France SAS - CE 2797



### INTELECT® MOBILE 2 COMBO:

| PART NUMBER   | DESCRIPTION   |
|---------------|---|
| 15-1204       | INTELECT® MOBILE 2 COMBO SET (EU PLUG)<br>(Combo set includes the below items)  |
| 15-1205       | INTELECT® MOBILE 2 COMBO SET (ALL PLUG)<br>(Combo set includes the below items) |
| 15-0133       | INTELECT® MOBILE 2 COMBO  |
| 79967         | Carbon electrodes   |
| 42198         | Gel electrodes  |
| 70010         | Stim lead wires   |
| 6522055       | Chattanooga® straps   |
| 15-0142       | 5cm <sup>2</sup> ultrasound applicator  |
| 4248          | Transmission gel bottle   |
| 15-0144/46/47 | Power cord  |

### INTELECT® MOBILE 2 STIM:

|               |   |
|---------------|---|
| 15-1202       | INTELECT® MOBILE 2 STIM SET (EU PLUG)<br>(Stim set includes the below items)  |
| 15-1203       | INTELECT® MOBILE 2 STIM SET (ALL PLUG)<br>(Stim set includes the below items) |
| 15-0132       | INTELECT® MOBILE 2 STIM   |
| 79967         | Carbon electrodes   |
| 42198         | Gel electrodes  |
| 70010         | Stim lead wires   |
| 6522055       | Chattanooga® straps   |
| 15-0144/46/47 | Power cord  |

### INTELECT® MOBILE 2 ULTRASOUND:

|               |   |
|---------------|---|
| 15-1200       | INTELECT® MOBILE 2 ULTRASOUND SET (EU PLUG)<br>(Ultrasound set includes the below items)  |
| 15-1201       | INTELECT® MOBILE 2 ULTRASOUND SET (ALL PLUG)<br>(Ultrasound set includes the below items) |
| 15-0131       | INTELECT® MOBILE 2 ULTRASOUND   |
| 15-0142       | 5cm <sup>2</sup> ultrasound applicator  |
| 4248          | Transmission gel (1 bottle)   |
| 15-0144/46/47 | Power cord  |

### INTELECT® MOBILE 2 OPTIONAL ACCESSORIES:

|         |   |
|---------|---|
| 15-1136 | Intelect Mobile 2 cart                  |
| 15-1210 | Intelect Mobile 2 cart with vacuum      |
| 70010   | Stim channel 1/2 lead wire kit          |
| 70012   | Stim channel 1/2 lead wire kit XL       |
| 15-1140 | USB drive                               |
| 79967   | 6x8cm carbon electrodes (4x)            |
| 15-0140 | 1cm <sup>2</sup> ultrasound applicator  |
| 15-0141 | 2cm <sup>2</sup> ultrasound applicator  |
| 15-0142 | 5cm <sup>2</sup> ultrasound applicator  |
| 15-0143 | 10cm <sup>2</sup> Ultrasound applicator |
| 4248    | Transmission gel (24 bottles)           |
| 14-1086 | Intelect Mobile 2 battery               |