

Intelect Mobile[®] 2 RPW USA

Radial Pressure Wave
technology made mobile



Intelect Mobile[®] 2 RPW



Chattanooga's Radial Pressure Wave technology is mobile.

The Intelect Mobile[®] 2 RPW (Radial Pressure Wave) is compact and light enough to be easily carried while helping to deliver a host of features and benefits.

When there is need for a powerful, yet portable RPW device in the clinic or on the sports field, the Intelect Mobile 2 RPW is an innovative choice. A lightweight and compact design together with a quiet built-in compressor make it suitable for portable use. A user-friendly experience is at the heart of the Mobile 2 RPW's design. The simple, ergonomically laid out controls enable all of the key parameters to be easily changed. A clear, illuminated LED display removes ambiguity when monitoring pressure, frequency and number of pulses. Pulses can be initiated via a trigger button on the handpiece.



User Profile

This Radial Pressure Wave clinical unit is ideally suited to the following applications:

- Orthopaedics
- Rheumatology
- Physical Medicine
- Physiotherapy
- Acupuncture
- Athletic trainers
- Podiatrists
- Chiropractors

Indications


- Reduce muscle pain and aches
- Temporarily increase blood flow
- Activate connective tissue

D-ACTOR^{®*} Applicator

The patented technology of the D-Actor[®] applicator is designed to aid in energy transmission at high frequency and intensive use. The therapist can easily connect a large choice of transmitters.



Standard D-ACTOR^{®*} Transmitters

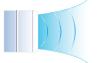
Ro40 	D 20-S 
0.38 mJ/mm² 	0.48 mJ/mm² 
0 - 35 mm	0 - 50 mm

Facts and figures

- Weight of the unit: Approx. 20 lbs
- Dimensions (L x W x H): 12.5 x 12.5 x 8"
- Built-in compressor
- Single and Continuous pulses
Frequency D-Actor: 1 - 21 Hz
Frequency V-Actor: 1 - 31 Hz
- Application pressure: 1 - 5 Bar
Highly customizable with precise adjustments of 0.1 increments
- Oscillating D-Actor®
- V-Actor handpiece option
- USB connectivity



User Interface Features

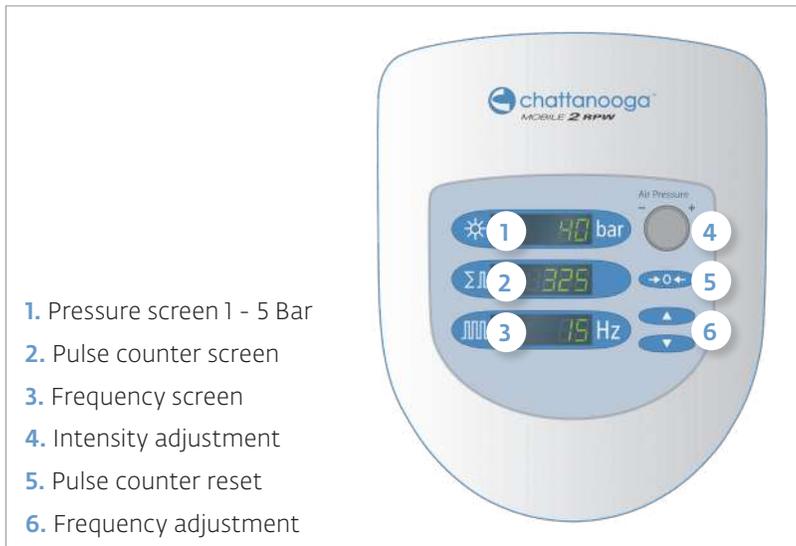

Hand piece docking port



Ultrasound gel holder



V-ACTOR®* Applicator

The V-ACTOR®* applicator is used for vibration massage therapy. When applied to tissue, the vibration pulse can help improve microcirculation, and enable elongation of fasciae and muscle fibers. It can also help to restore normal muscle tone.



Standard V-ACTOR®* Transmitters

