


TOM 4 Classic

Fully adjustable and foldable budget rehab-stroller



- durable and approved
- full positioning features
- stabilization in sitting for everyday activities.

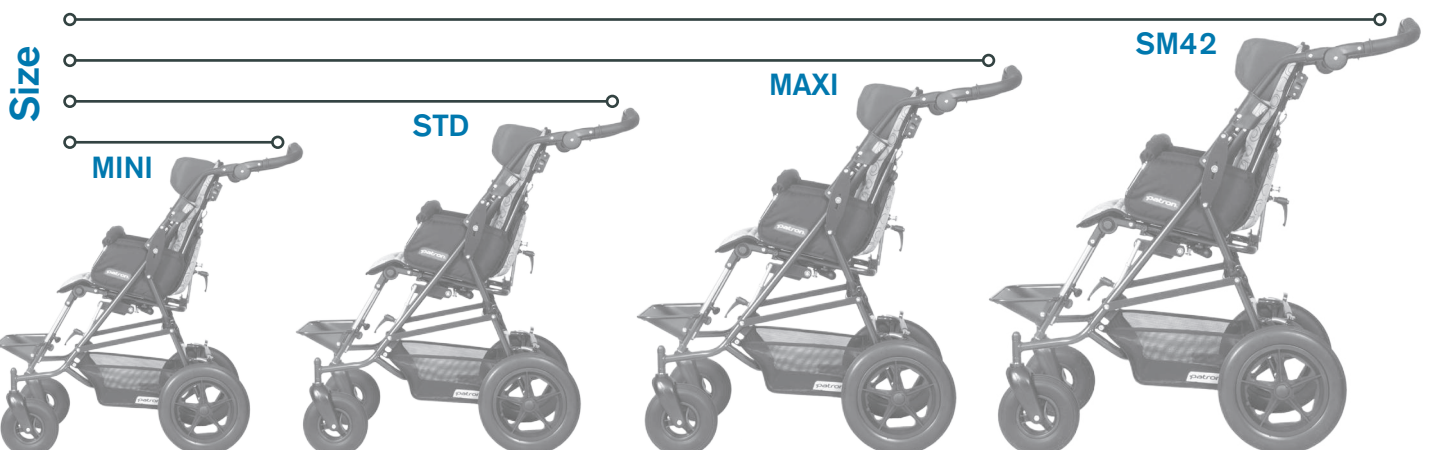

TOM 4 Classic – product line of fully adjustable rehab-strollers for handicapped children is based on improved design of the TOM-line undercarriages with 200/300 mm solid wheels, that enables usage of the product in town and as well as for the outdoor activities.

Fully foldable undercarriage is produced from lightweight thin-walled tubes with powder coating finish. Redesigned seat-unit is detachable, reversible and fully size and position adjustable. It includes as well new flexible system for attaching of fixation devices and accessories. Seat shell is covered by fully detachable soft padded textile upholstery in fashionable design. Easy handling is ensured as well as by the stroller's fully foldable frame even with attached seat. All the products correspond to standard hygienic, technical and safety requirements.

Standard equipment:

- detachable and reversible seat-unit
- height adjustable backrest and depth adjustment of the seat-base.
- width adjustable seat by FlexiPads set (optional equipment)
- seat with fully adjustable tilt in space function
- fully reclining backrest
- angle and height adjustable footrest
- 200/300 mmsolid wheels
- detachable soft side-cushions
- detachable grab bar with cover
- frame colour and fabric patterns according to the current offer



folded stroller

Optional equipment:



Canopy CLASSIC



Basket (textile)



Work board BASIC - UNI



FLEXI-pads set



Rain cover



Headrest CLASSIC



Side trunkrest CLASSIC



4-point immobilisation belt CLASSIC



Hip rest CLASSIC



Abduction block CLASSIC



Immo. Waistcoat CLASSIC



H-style harness CLASSIC



Pelvis fixation



Feet fixation - pair

.... and
much more
accessories!

according to the actual list

Full scale of growing functions




Full seat positioning



fully adjustable



reversible seat

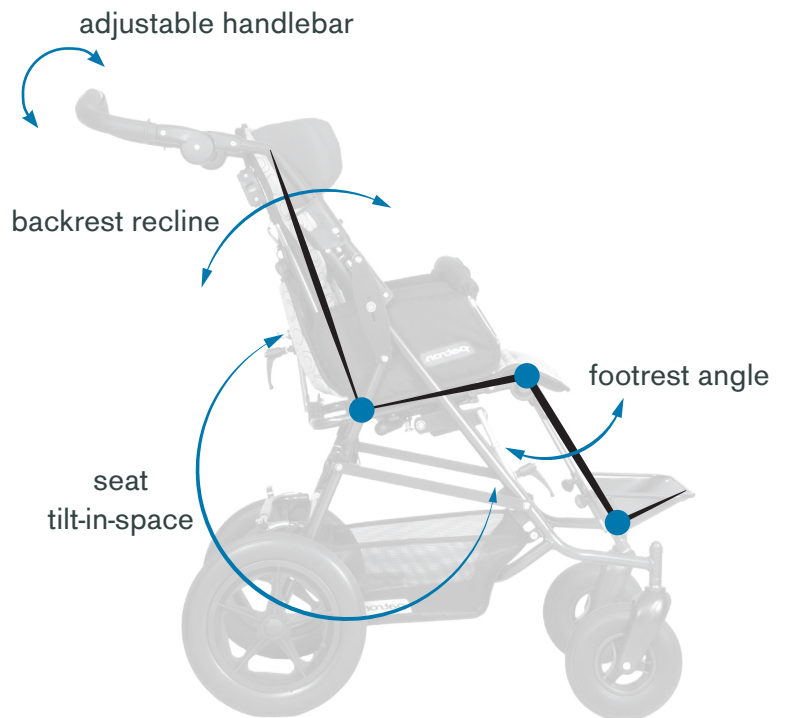

DESIGN OPTIONS



R43



R44



R45



R50

SUCCESSFULLY CRASH TESTED

The product had been successfully crash tested and fulfilled all the requirements of EN ISO 7176/19. It enables to use the rehab-buggy not only for a common duty, but as well as for the transportation in specially adjusted motor vehicles, equipped with the anchoring points and equipment.



Order codes:

TOM 4 Classic MINI:	T4CWKPMYY
TOM 4 Classic STD:	T4CWKPYYY
TOM 4 Classic MAXI:	T4CWKPXY Y
TOM 4 Classic SuperMAXI 42 (reinforced version):	T4CWKPS2Y

	MINI	STD	MAXI	SuperMAXI 42
Backrest height [mm]	450-600	520-680	570-680	570-680
Seat depth [mm]	270-310	270-310	320-430	320-430
Seat width [mm]	260/300	240/280/340	280/320/380	425
Seat to footrest [mm]	140-280	200-350	240-400	240-400
Backrest angle	0° - 175°	0° - 175°	0° - 165°	0° - 165°
Seat angle	+/-27°	+/-27°	+/-25°	+/-25°
Total length [mm]	990	990	1070	1100
Total width [mm]	580	580	630	680
Total height /with adj.handle [mm]	1010	1010	1060	1030
Brake	foot oper.	foot oper.	foot oper.	foot oper.
Weight	17,0 kg	18,0 kg	20,5 kg	24,0 kg
Maximum load	35,0 kg	35,0 kg	50,0 kg	75,0 kg
Folded stroller length [mm]	930	930	960	960
Folded stroller width [mm]	580	580	630	680
Folded stroller height [mm]	550	550	630	660